

SURMODICS VASCULAR DEVICE SOLUTIONS

Low Profile .035" Balloon Dilation Catheters

.035"
LOW PROFILE

ULTRA LOW PROFILE FOR EXCEPTIONAL DELIVERABILITY

DESIGNED TO NAVIGATE
THE MOST CHALLENGING
VASCULAR ANATOMY




SURMODICS

Unmatched low profile


- 5Fr sheath compatibility through 8mm balloon diameters
- State-of-the-art balloon pleat & fold process ensures superb crossing profiles

Exceptional delivery in challenging vascular anatomy

- Hydrophilically coated for smooth tracking and crossability
- Dual lumen shaft designed for kink resistance and responsive push
- Laser polished tip provides gradual transitions for lesion entry

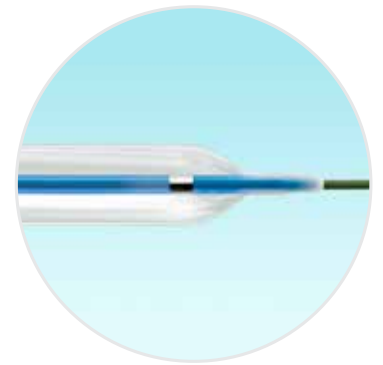

Product Features

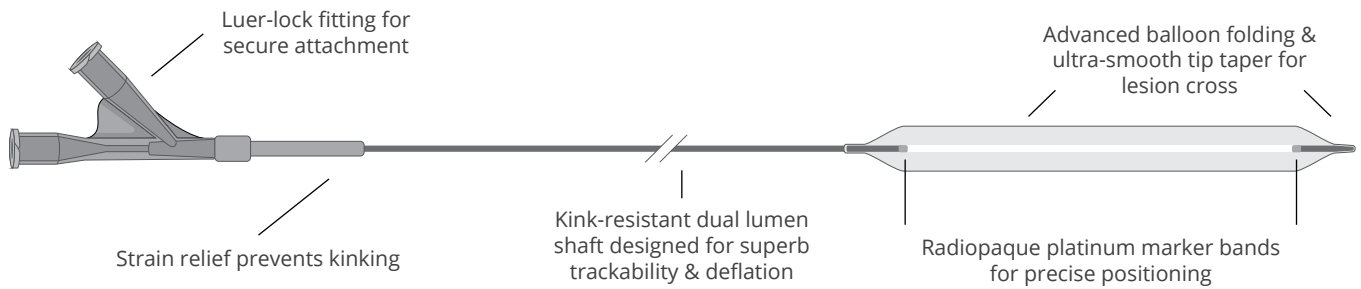

Product Specifications

Catheter Platform	OTW
Recommended Guide Wire	.035"
Sheath Compatibility	5Fr (through 8mm x 60mm)
Catheter Working Lengths	80cm & 135cm (9mm x 80mm & 10mm x 80mm is 80cm only)
Balloon Diameters	5mm - 10mm
Balloon Lengths	20mm - 150mm (9mm & 10mm are 20mm - 80mm)
Nominal Pressure	8atm
Rated Burst Pressure	12atm (5mm - 8mm) 10atm (9mm & 10mm)



Creagh Medical and NorMedix are now part of **Surmodics, Inc.**

9924 West 74th Street
Eden Prairie, MN 55344 USA
Toll Free 1-866-787-6639
Phone 952-500-7000
Fax 952-500-7001
www.surmodics.com

© 2016 Surmodics, Inc.
All rights reserved. 35LP0416 01

FOR MORE INFORMATION ON HOW **SURMODICS** CAN ENHANCE YOUR BUSINESS, PLEASE CONTACT US AT **952-500-7000** OR VISIT **WWW.SURMODICS.COM**.